












Biocard®

Malaria Pf/Pv Rapid Test Kit

Biocard® Malaria Pf/Pv Rapid test is an immunochromatographic assay for the qualitative detection of *Plasmodium falciparum* (Pf) and *Plasmodium vivax* (Pv) antigen in the human whole blood specimen. The test band (Pf) is pre-coated with anti HRP2 monoclonal antibody and the test band (Pv) is pre-coated with anti Pv pLDH monoclonal antibody and a control line with control antibody.

-  Technique: Rapid test based on lateral flow
-  Accuracy: Meets WHO recommended detection criteria
-  Detection: *Plasmodium vivax* and *Plasmodium falciparum*
-  Results: 20 minutes
-  Longer shelf life: 24 months
-  Storage: 4°C to 30°C
-  Convenient pack size: 30 tests
-  Sensitivity: 100%
-  Specificity: 100%


Easy to Use

- Simple test using only 1 buffer reagent
- Ease of use, requires minimal training of personnel
- Requires only 5 µL of venous or capillary whole blood
- No high cost equipment required for testing
- No need for skilled microscopist


Easy to Interpret

- Meets WHO acceptance criteria
- Excellent correlation with microscopic results
- Differentiated results allow diagnosis of potentially lethal *P. falciparum* infection
- Easy to read, qualitative results
- Higher Specificity and better Sensitivity

Test Procedure

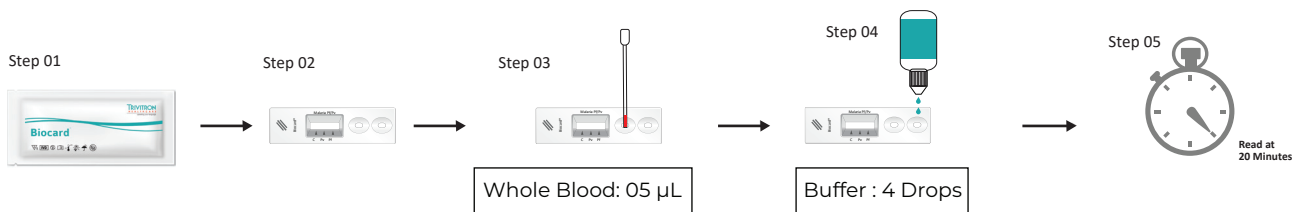

Result Interpretation

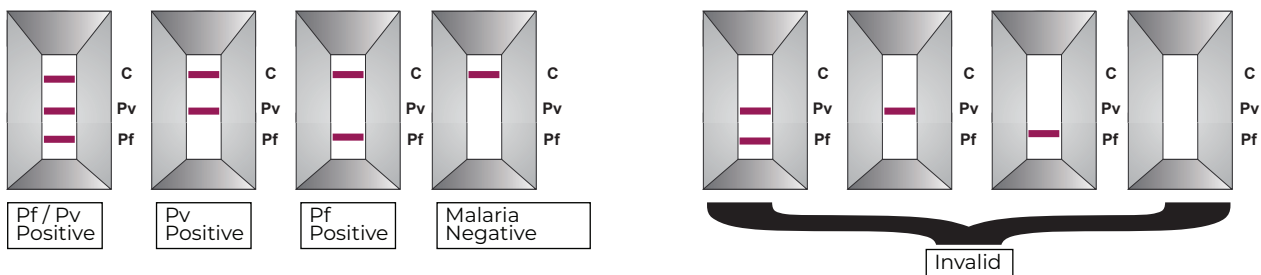

Order Information

Catalog	Product	Pack Size	Specimen
30201437	Biocard® Malaria Pf/Pv Rapid Test	30 Tests	Whole Blood

*Please note: Product specifications are subject to change without prior notice owing to product modifications, improvements / up-gradation.



corporate@triviron.com | www.triviron.com

Customer Care Helpline: +91 98400 80008



TRIVITRON
HEALTHCARE
speaking your language